

SIGER

**SIGER
EMPOWERING FUTURE**

V6000 NEW
Dental Unit

THE BEAUTY
OF SIMPLICITY



Zhuhai Siger Medical Equipment Co., Ltd

Add : Building 2, No.1 Chuangxin Yi Road, Tangjiawan Town, Zhuhai City, Guangdong Province, P. R. China

Tel : +86 756 3881087 Email : info@siger.cn Web : en.siger.cn

V6000

Dental Unit

Simple yet extraordinary

- New design of delivery system, **Innovative spatial layout**
- Brand new **intelligent control system** with 3 doctor programs, making personalized settings possible
- Compact packaging design permits **efficient installation**
- **Easily fits the entire package** into the clinic
- Ergonomically designed dental chair with ultra-soft upholstery, **provides an excellent comfort experience**



STRUCTURAL CONFIGURATION

A Brand New Delivery System

Specially designed for higher efficiency and comfort

The treatment table can be flexibly positioned at a wide angle

Meets diverse needs of working position

● Operation activity zone

NEW LAYOUT

Redefined Layout and Space

Striving for excellent treatment experience

Reasonable spatial zoning for doctors and patients to improve efficiency

Ultimate comfort, relaxed treatment

● Patient activity zone

TREATMENT TABLE

Intelligent Control System

Base-mounted treatment table,
flexible moving, and stable positioning

Double joint & three-axis delivery arm



- Delivery arm directly connected to the chair base plate, combining agility and stability.
- Optional height-adjustable arm, catering to multiple work demands.

Powerful and concise doctor control panel



- Touch buttons, clear icons, straightforward operation.
- 3 doctor programs + 15 chair memory positions, allowing multi-user personalization.
- Smart linkage program for cup filling, spittoon rinsing, chair position and dental light.
- Thoughtful safety settings, one-button stop function and interlock protection.

Extra large treatment table with 6 instrument slots



- Large space, strong bearing support.
- 6 instruments available, compatible with various dynamic instruments.



Newly Designed Dental Chair Integration of Comfort and Esthetics

- Upgraded chair structure brings smooth and solid movements.
- Unique SmoothControl™ chair comforting system, makes every movement stable.



Ergonomic design for enhanced comfort

- Well-made upholstery, provides excellent support from head to legs.
- Evenly distributes the weight of the body, helps to keep relaxed during long treatment.



Combination of safety, practicality, and comfort

- Ultra-thin backrest, creates larger legroom for doctor and assistant.
- Low minimum position, convenient for children and elderly patients.
- Angle sensor technology used on chair position memory, ensuring accuracy, avoiding cumulative error.
- Chair descent protection to prevent accidents.




Premium materials and advanced manufacturing process

- Standardized production process for molds: plastic injection and free from wood materials, promising hygiene and durability.
- Snap-fit upholstery system, customizable colors to perfectly match clinic design.
- High-strength steel, ensuring stability and outstanding rust prevention performance.

MULTIFUNCTIONAL FOOT CONTROL

Pedal/lever type available

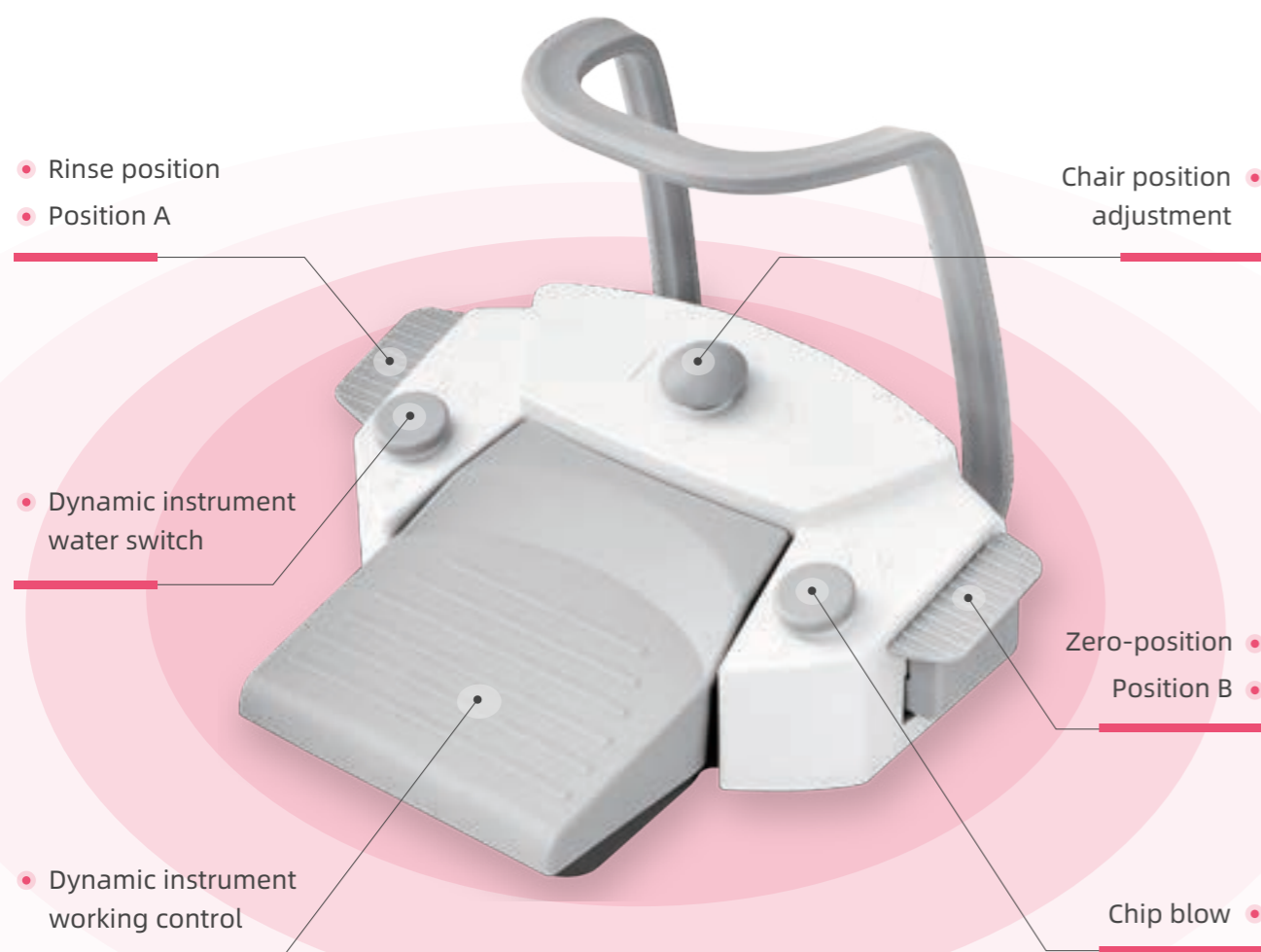

Various Options on Dental LED Lights

Blue light cutoff
Adjustable illuminance



SIGER V1

Light source: LED
Color temperature: 4300K
Illuminance: 8000-30000LUX
Light spot: 70x140MM(@700MM)
Power: 15VA
Voltage: 12-24Vac



SIGER V2

Light source: LED
Color temperature: 5000K
Illuminance: 8000-30000LUX
Light spot: 80x150MM(@700MM)
Power: 6VA
Voltage: 12-24Vac



SIGER V5

Light source: LED
Color temperature: 3000-5500K
Illuminance: 8000-35000LUX
Light spot: 90x200MM(@700MM)
Power: 18VA
Voltage: 12-24Vac




L500

Light source: LED
Color temperature: 3300-4300K
Illuminance: 33000-130000LUX
Light spot: 120x200MM(@1000MM)
Power: 25VA
Voltage: 12-24Vac

NEW PACKAGING DESIGN

Ingenious design of product and packaging

Unique packaging and product design for easy installation and maintenance

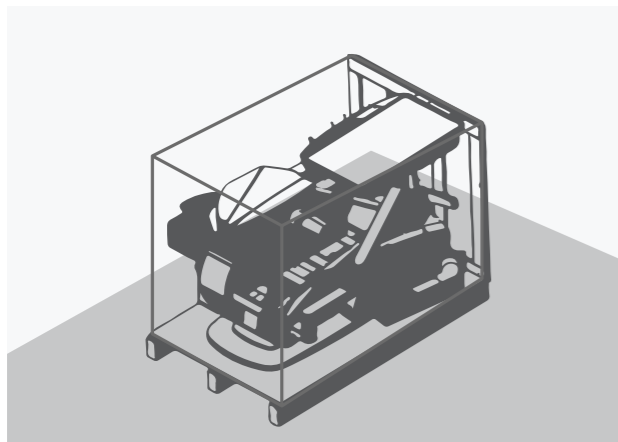

- ▶ Compact packaging box, ensuring transportation safety while greatly enhancing entry convenience.



- ▶ Carefully designed components and structure, make installation easier.



- ▶ Mainly adopted plug-in structure from part to part, minimal tool usage, greatly improving efficiency.



- ▶ Efficient utilization of space, saving logistics and transportation costs.

EXPRESS SETUP

Quick Installation

Minimalist packaging for time-saving and effortless installation.



EXTRAORDINARY QUALITY

Outstanding quality & careful consideration of hygiene

Thoughtful design with high quality



- Detachable suction tubes and filters for daily cleaning



- Detachable and disinfectable spittoon bowl and water outlet



- Adding distilled water from the top: Easy and convenient



- ABS plastics, standardized molding process



- High-quality and durable pipelines



- Excellent painting, good rust-proof performance

OTHER CONFIGURATION

Stylish and colorful

Standard foam upholstery



Optional sewn upholstery (silicone-coated leather)



V6000 dentist stool

- ▶ The seat height and the backrest's tilt can be adjusted to meet ergonomics.
- ▶ Noise-reducing casters.



• Optional assistant stool

THE ACHIEVEMENT OF ERA CRAFTSMANSHIP

Crystallization of modern aesthetic connotation

User-friendly design, serve for doctor and patient

